



Leica EM GP

Automatic Plunge Freezer for the Bare Grid Technique

Living up to Life

Leica
MICROSYSTEMS

The Bare Grid Technique

Many specimens for cryo-TEM can be prepared by immersion freezing, where a fluid sample is pipetted onto an EM grid (usually coated) and the excess removed until a thin film remains, before plunging into a cryogen such as liquid ethane. The grid can then be directly transferred under cryo conditions to the cryo electron microscope (cryo-TEM) for observation. This is the bare grid technique.

The bare grid technique can be used for many types of sample ranging from biological cell sub units to industrial emulsions. Imaging macromolecular assemblies, viruses and cells in their native, hydrated environment in the cryo-TEM is the state-of-the-art technique in electron microscopy, providing maximum resolution with minimal specimen damage.

Although a simple method, it is imperative that the suspension thickness on the grid can be reproduced and vitreous ice can be formed, otherwise much time is wasted loading useless samples into the cryo-TEM. The sample film is only tens to hundreds of nanometers thick and so can be easily influenced by temperature shifts and humidity prior to freezing. If the humidity is too low then the film breaks due to it quickly drying. If the sample is adversely affected by temperature, then the morphology may change before freezing.

Leica Microsystems has developed a plunge freezer, in conjunction with Dr. Guenter Resch of the IMP/IMBA Electron Microscopy Facility in Vienna, Austria, to standardize procedures and make the bare grid technique more reproducible.



Leica EM GP Form and Function



The Leica EM GP plunge freezes samples into a secondary cryogen such as liquid ethane after removing excess fluid by automatic blotting.

After connecting the forceps holding the grid to the Leica EM GP, an environmental chamber lowers to surround the grid, providing a protective temperature and humidity controlled environment.

Access ports on both sides of the environmental chamber allow easy pipetting of solutions and suspensions for both left and right handed users.

Excess fluid is then removed by automatic blotting with filter paper from one side of the grid. Parallel, single sided blotting was developed to prevent damage to delicate support films. The blotter touches the complete grid surface in one movement. The grid can be programmed to automatically turn 180° prior to blotting to allow fluid removal from the correct side.

All parameters are displayed on the touch-screen control panel, where numerous settings can be adjusted and monitored such as blot time, grid/blotter positioning, temperature, humidity, LN₂ level, and secondary cryogen temperature.

After freezing, the grid is transferred to a pre-cooled grid box inside a transfer container filled with LN₂ located in the Dewar. This can then be taken to the sample holder of the cryo-TEM to prepare for loading into the electron microscope.

Leica Design by Werner Hölbl

Leica EM GP Features

Environmental chamber

To provide the best conditions for pre-freezing of a suspension, the environment surrounding the grid and sample has to be precisely controlled. The environmental chamber envelops the grid and forceps, providing a temperature and humidity controlled protective environment variable between +4° and +60°C and room humidity to 99%.


To provide a clear view, an anti-fogging heater keeps the glass window clear.

Blotting

Excess fluid is removed by automatic blotting which can be initiated in two ways:

- 1. Aut blotting** – after the suspension has been pipetted onto the grid and the button pressed, the blotter automatically moves to the pre-set blot position for the user-programmed time before plunging of the grid.
- 2. Sensor control** – the most automated method of plunging. After pressing the plunge button the blotter moves towards the grid. A photosensor detects the moment the droplet touches the grid and the blotting time countdown begins. The grid is then either held for the desired length of time (to allow redistribution of fluid for example) or immediately plunged into the secondary cryogen.

Blotting and plunging can be activated via the touch screen or footswitch.

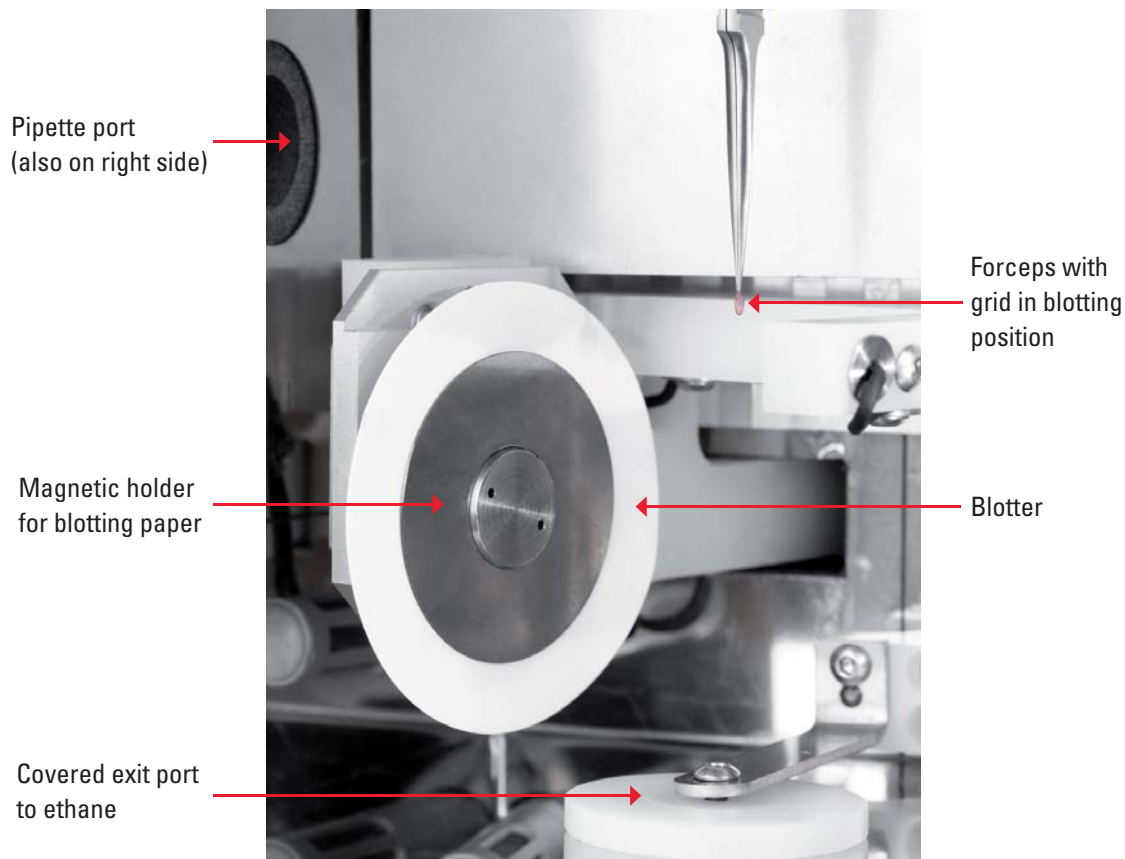

Viewing system

The Leica EM GP has an optional stereomicroscope to aid alignment and sample preparation. The LED illumination, fixed both inside the environmental chamber and beaming down onto the Dewar, provides excellent high light levels for observation of the complete preparation and plunging process.

The Dewar

After switching on the Leica EM GP, the 1 liter Dewar can be filled with LN₂ before liquefying the secondary cryogen, usually ethane. A full Dewar lasts for approximately 1 hour between refilling. Liquefying the secondary cryogen is fast, easy and safe with the liquefying head. The head is connected to the secondary cryogen regulator of the gas bottle and the gas slowly fed in. It condenses within seconds, taking about one minute to fill the 2.5 ml container. A cover is provided to prevent LN₂ splashing into the ethane on subsequent refilling of the LN₂. The temperature of the secondary cryogen can be controlled precisely from the control panel. A container, filled with LN₂, sits in the Dewar to hold a grid box for transfer of prepared, vitrified samples.



Liquefier in place over ethane container in Dewar



After freezing the grid remains in or above the ethane (depending upon user settings) ready for transfer to the grid box

Control Panel

Operation is via touch screen control with all adjustable parameters visible.

Program screen

All parameters can be adjusted and set for up to 10 programs. In the Setup menu the positioning of the grid relative to the blotter can be adjusted and also the transfer position after plunging.

Safety

The Leica EM GP operates under strict safety conditions. During any movement of the environmental chamber a large red STOP

button appears on the control panel. Touching this button will immediately stop any movement. An alarm signals when either the secondary cryogen is too warm and may evaporate or the LN₂ level is too low.

Bake-out

At the end of a run the bake-out cycle takes 60 minutes to dry the Dewar and environmental chamber, which allows a second run within a short time if the user does not wish to maintain the LN₂ level in the Dewar.

Main screen

Program in use → P 2 5.0/ .5/ 1.0/

Current temperature of environmental chamber → T_c 37°C

Current humidity of environmental chamber → H_r 80 %

Current temperature of ethane container → LN₂ -180°C

LN₂ level → 50%

Count → 10

Light → [Sun icon]

STOP

Load Forceps → Sets start position to accept forceps

Lower Chamber → Chamber moves downwards ready for dispensing suspension

Rotate Home Blot /A-Plunge → Blotting/plunging sequence

Plunge

Transfer → Press to raise grid in ethane and increase Dewar GN₂ production to keep out moisture

Setup Bake out Manual Rotation

Countdown for blotter exchange

Program screen

Nr.	Blotting Sensor	A-Plunge	Delay	Blot	Hold
1	✓	✓	0.0	0.0	0.0
2	✓	✓	0.0	0.5	0.0
3	✓	✓	0.0	0.0	0.0
4	✓	✓	0.0	0.0	0.0
5	✓	✓	0.0	0.0	0.0

Rotate grid after applying suspension

Plunge automatically after blotting

Blotting time

Program number → 1

Scroll through programs 1-5 or 6-10

Rotate grid before applying suspension

Blotting with blot sensor activated

Delay time before blotting


Delay before freezing after blotting

OK


Applications

The Leica EM GP is designed for all EM laboratories with a need to view vitrified fluid samples or extremely thin samples in the cryo-TEM, including biological research, virology, protein crystallography, pharmaceutical research, cosmetics and industrial laboratories.

Samples that can be prepared vary for example from suspensions of viruses, liposomes, microtubules, proteins and other cellular components to paint or solutions and emulsions in both aqueous and inorganic solvents. The Leica EM GP can be used to plunge freeze samples not only on EM grids for the Bare Grid Technique, but also sapphire discs and samples in freeze fracture planchettes.

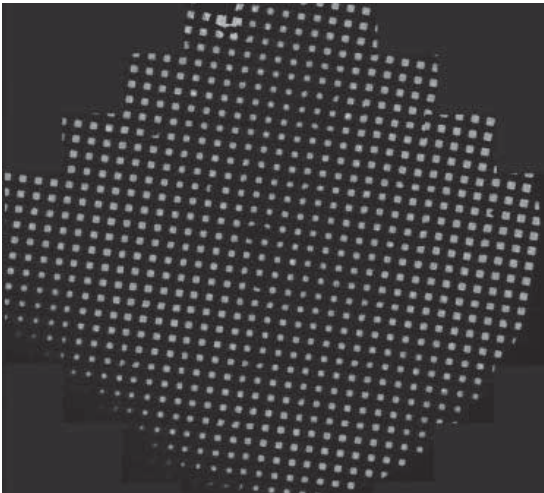


Montage overview of plunge frozen grid.
Note the homogeneity of the film thickness.

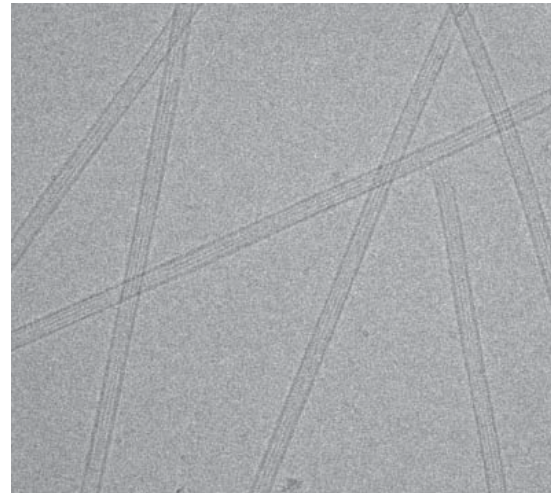


Microtubules

Micrographs courtesy of Dr. Guenter Resch, IMP/IMBA Electron Microscopy Facility, Vienna, Austria

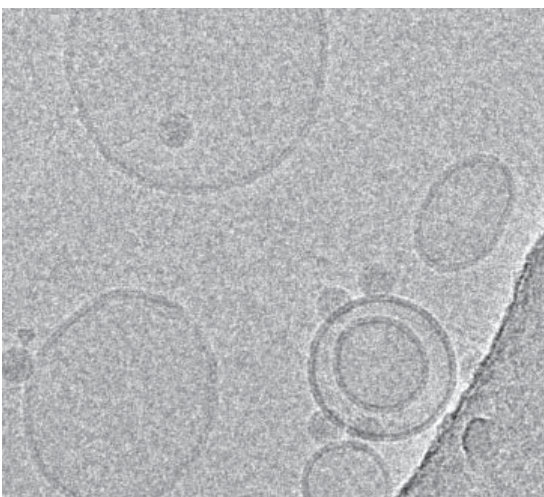

Liposomes

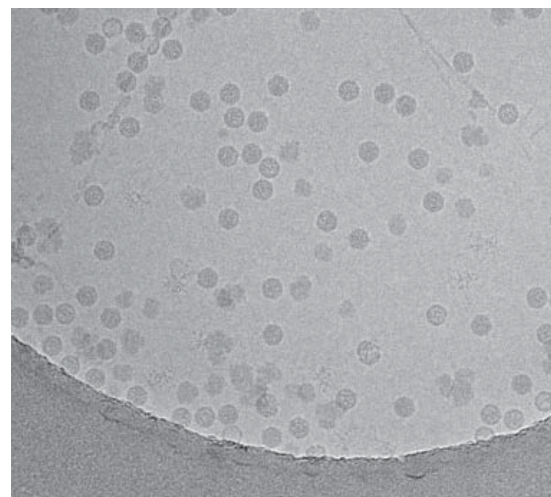

Rhinovirus particles on holey carbon film

Micrographs courtesy of Angela Pickl-Herk, MFPL, Vienna, Austria

“With the user, for the user”

Leica Microsystems

Leica Microsystems operates globally in four divisions, where we rank with the market leaders.

• Life Science Division

The Leica Microsystems Life Science Division supports the imaging needs of the scientific community with advanced innovation and technical expertise for the visualization, measurement, and analysis of microstructures. Our strong focus on understanding scientific applications puts Leica Microsystems' customers at the leading edge of science.

• Industry Division

The Leica Microsystems Industry Division's focus is to support customers' pursuit of the highest quality end result. Leica Microsystems provide the best and most innovative imaging systems to see, measure, and analyze the microstructures in routine and research industrial applications, materials science, quality control, forensic science investigation, and educational applications.

• Biosystems Division

The Leica Microsystems Biosystems Division brings histopathology labs and researchers the highest-quality, most comprehensive product range. From patient to pathologist, the range includes the ideal product for each histology step and high-productivity workflow solutions for the entire lab. With complete histology systems featuring innovative automation and Novocastra™ reagents, Leica Microsystems creates better patient care through rapid turnaround, diagnostic confidence, and close customer collaboration.

• Surgical Division

The Leica Microsystems Surgical Division's focus is to partner with and support surgeons and their care of patients with the highest-quality, most innovative surgical microscope technology today and into the future.

The statement by Ernst Leitz in 1907, “with the user, for the user,” describes the fruitful collaboration with end users and driving force of innovation at Leica Microsystems. We have developed five brand values to live up to this tradition: Pioneering, High-end Quality, Team Spirit, Dedication to Science, and Continuous Improvement. For us, living up to these values means: **Living up to Life.**

Active worldwide

Australia:	North Ryde	Tel. +61 2 8870 3500	Fax +61 2 9878 1055
Austria:	Vienna	Tel. +43 1 486 80 50 0	Fax +43 1 486 80 50 30
Belgium:	Groot Bijgaarden	Tel. +32 2 790 98 50	Fax +32 2 790 98 68
Canada:	Richmond Hill/Ontario	Tel. +1 905 762 2000	Fax +1 905 762 8937
Denmark:	Herlev	Tel. +45 4454 0101	Fax +45 4454 0111
France:	Nanterre Cedex	Tel. +33 811 000 664	Fax +33 1 56 05 23 23
Germany:	Wetzlar	Tel. +49 64 41 29 40 00	Fax +49 64 41 29 41 55
Italy:	Milan	Tel. +39 02 574 861	Fax +39 02 574 03392
Japan:	Tokyo	Tel. +81 3 5421 2800	Fax +81 3 5421 2896
Korea:	Seoul	Tel. +82 2 514 65 43	Fax +82 2 514 65 48
Netherlands:	Rijswijk	Tel. +31 70 4132 100	Fax +31 70 4132 109
People's Rep. of China:	Hong Kong	Tel. +852 2564 6699	Fax +852 2564 4163
Portugal:	Lisbon	Tel. +351 21 388 9112	Fax +351 21 385 4668
Singapore		Tel. +65 6779 7823	Fax +65 6773 0628
Spain:	Barcelona	Tel. +34 93 494 95 30	Fax +34 93 494 95 32
Sweden:	Kista	Tel. +46 8 625 45 45	Fax +46 8 625 45 10
Switzerland:	Heerbrugg	Tel. +41 71 726 34 34	Fax +41 71 726 34 44
United Kingdom:	Milton Keynes	Tel. +44 1908 246 246	Fax +44 1908 609 992
USA:	Bannockburn/Illinois	Tel. +1 847 405 0123	Fax +1 847 405 0164

and representatives in more than 100 countries